

CAJA INTERSOMÁTICA



Caja intersomática (TLIF)
Interbody Cage (TLIF)

Características Técnicas

Technical specifications

Las Cajas de Surgival permiten una artrodesis intersomática por vía transforaminal unilateral con buenos resultados clínicos y con escasa incidencia de complicaciones.

Longitudes:
26, 30 mm

Alturas:
7, 8, 9, 10, 12 y 14 mm

Material:

- Marcadores radio-opacos de Tantalio (R05400).
- Fabricado en Peek Optima®.

Morfología:

- Superficie cráneo-caudal ranurada con ángulo de lordosis de 4°.
- La geometría en forma de luna creciente aumenta la estabilidad.
- El diseño permite la transmisión de cargas al injerto óseo.
- El nervio central biconcavo aumenta la solidez de la caja ante fuerzas de cizallamiento sin disminuir la superficie de exposición del injerto óseo que contiene en su interior.

Vía de abordaje transforaminal (TLIF):

- Debe asociarse siempre con una fijación posterior.
- Sólo debe utilizarse un implante por nivel a artrodesar.

Surgival Cages are used to perform lumbar interbody arthrodesis using a unilateral transforaminal with good clinical results and scarce incidents of complications.

Lengths:
26, 30 mm

Heights:
7, 8, 9, 10, 12 and 14 mm

Material:


- Tantalum (R05400) radio-opaque markers to control.
- Manufactured in Peek Optima®.

Morphology:

- Grooved cranio-caudal surface with 4° lordosis angle.
- The geometry in the shape of a crescent moon increases stability.
- The design enables the transmission of loads on bone graft.
- The central bioncave nerve increases the solidity of the cage against shearing forces without diminishing the bone graft's exposure surface, which is contained inside it.

Transforaminal Lumbar Interbody Fusion (TLIF):

- It must always be associated with posterior fixation.
- Only one implant per level to be arthrodesed must be used.

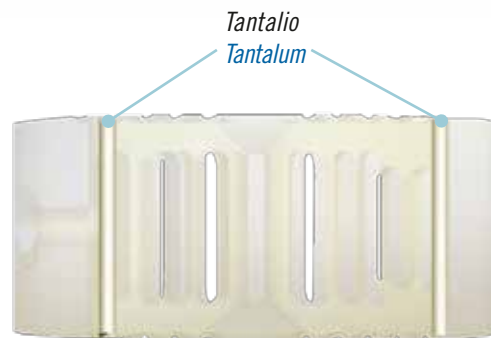



Caja intersomática (TLIF) Interbody Cage (TLIF)		
	Longitud Length	Altura Height
Ref. C6890007E	26 mm	7 mm
Ref. C6890008E	26 mm	8 mm
Ref. C6890009E	26 mm	9 mm
Ref. C6890010E	26 mm	10 mm
Ref. C6890012E	26 mm	12 mm
Ref. C6890014E	26 mm	14 mm
Ref. C6890207E	30 mm	7 mm
Ref. C6890208E	30 mm	8 mm
Ref. C6890209E	30 mm	9 mm
Ref. C6890210E	30 mm	10 mm
Ref. C6890212E	30 mm	12 mm
Ref. C6890214E	30 mm	14 mm

CAJA INTERSOMÁTICA INTERBODY CAGE

Set completo instrumental Caja intersomática
Interbody Cage instrument complete set
Ref. C6890000



Set 1 instrumental Caja Intersomática
Interbody Cage instrument set 1

1 Distractor
Distractor
Ref. C6892435

2 Cureta recta
Straight curette with serrated cup
Ref. C6892375
Cureta curvada derecha
Right curved curette with serrated cup
Ref. C6892380
Cureta curvada izquierda
Left curved curette with serrated cup
Ref. C6892385



3 Legra recta
Straight curette with hollow cup
Ref. C6892360
Legra curvada derecha
Right curved curette with hollow cup
Ref. C6892365
Legra curvada izquierda
Left curved curette with hollow cup
Ref. C6892370

4 Escoplo Chisel	Longitud Length
Ref. C6892308	8 mm
Ref. C6892312	12 mm

Instrumental Instruments

Set 2 instrumental Caja Intersomática Interbody Cage instrument set 2

1	Escoplo cuadrangular Box chisel	Altura Height
	Ref. C6892207	7 mm
	Ref. C6892208	8 mm
	Ref. C6892209	9 mm
	Ref. C6892210	10 mm
	Ref. C6892212	12 mm
	Ref. C6892214	14 mm

2	Raspa de prueba 26 mm Trial rasp 26 mm	Altura Height
	Ref. C6891007	7 mm
	Ref. C6891008	8 mm
	Ref. C6891009	9 mm
	Ref. C6891010	10 mm
	Ref. C6891012	12 mm
	Ref. C6891014	14 mm

	Raspa de prueba 30 mm Trial rasp 30 mm	Altura Height
	Ref. C6891207	7 mm
	Ref. C6891208	8 mm
	Ref. C6891209	9 mm
	Ref. C6891210	10 mm
	Ref. C6891212	12 mm
	Ref. C6891214	14 mm

3	Separador para saco dural Separator for dural sac
	Ref. C6892350
	Separador para raíz nerviosa Separator for nerve root
	Ref. C6892355

6	Casquillo para fresas Bushing for reamers
	Ref. C6892450



8	Fresa manual Manual reamer	Altura Height
	Ref. C6892107	7 mm
	Ref. C6892108	8 mm
	Ref. C6892109	9 mm
	Ref. C6892110	10 mm
	Ref. C6892112	12 mm
	Ref. C6892114	13 mm

4	Mango intercambiable en T para fresas Interchangeable T handle for reamers
	Ref. C6892010

7	Barra extractora para mango intercambiable de anclaje rápido Extractor bar for quick coupling interchangeable handle
	Ref. C6892006

5	Mango intercambiable anclaje rápido Quick coupling interchangeable handle
	Ref. C6892005


Set 3 instrumental Caja Intersomática Interbody Cage instrument set 3

1	Varilla contención de Caja Cage contention rod
	Ref. C6892416

2	Introducir de Caja Cage inserter	Longitud Length
	Ref. C6892410	26 mm
	Ref. C6892415	30 mm

3	Empujador recto Straight pusher	Altura Height
	Ref. C6892440	6 mm
	Ref. C6892445	10 mm
	Empujador curvado Curved pusher	Altura Height
	Ref. C6892420	6 mm
	Ref. C6892425	10 mm

4	Rescatador de Caja Cage salvor
	Ref. C6892430

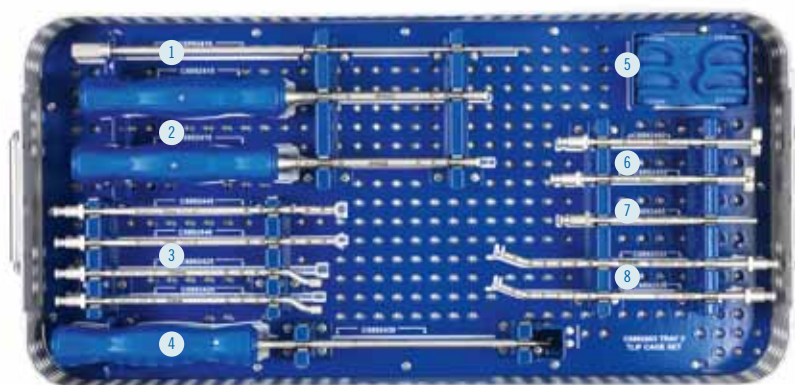

5	Soporte para Caja Cage support
	Ref. C6892400

6	Compactor grande de injerto Large graft compactor	Longitud Length
	Ref. C6892455	26 mm
	Ref. C6892460	30 mm

7	Compactor pequeño de injerto Small graft compactor
	Ref. C6892405

8	Empujador angulado Angled pusher	Altura Height
	Ref. C6892530	6 mm
	Ref. C6892535	10 mm



Tecnología Europea de Vanguardia
Advanced European Technology

Distribuido por:

Distributed by:



SURGIVAL COMERCIAL

C/Ignacio Iglesias, 70
08950 Esplugues de Llobregat · Barcelona · España
Tel: (+34) 93 480 92 22 · Fax: (+34) 93 480 92 23
e-mail: orbimed@surgival.com

FÁBRICA / DEPARTAMENTO INTERNACIONAL

Parque Tecnológico
C/ Leonardo Da Vinci, 12-14 · 46980 Paterna · Valencia · España
Tel: (+34) 96 131 80 50 · Fax: (+34) 96 131 80 95
e-mail: surgival@surgival.com

www.surgival.com

